

DeNovix®

CellDrop™

Count Cells Without Slides™



SCIENTISTS' CHOICE AWARDS 20
WINNER

Patented
US 10,302,568
US 10,598,598
JP 6579475

Independent reviews from SelectScience.net

“ Game changing. ”

– A. Rousseau, MRC PPU

Rating: 5.0 ★★★★★

“ Every cell room should have one! ”

– F. Fomeris, University of Pavia

Rating: 5.0 ★★★★★

“ CellDrop is a great cell counter! ”

– T. Smith, Duke University

Rating: 5.0 ★★★★★



Reviewers' Choice
Life Sciences Customer Service
Company of the Year

Brightfield and Fluorescence Automated Cell Counters

Count Cells Without Slides™

DeNovix CellDrop™ Automated Cell Counters combine patented DirectPipette™ sample loading with advanced fluorescence and brightfield optics as the new standard in cell counting.

The CellDrop Family removes the need for disposable slides from routine cell counting bringing low-volume, consumable-free analysis to the cell biology lab.

CellDrop received the Best New Life Science Product of the Year at the 2020 Scientists' Choice Awards®.




EasyApps® software comes preinstalled with cell counting and viability apps for AO/PI, Brightfield, Trypan Blue, GFP, Yeast as well as powerful tools for creating custom applications.

Eliminate Slide Costs


- No Purchase Costs
- No Slide Disposal Fees
- Never Run Out of Stock
- No Shipping or Storage Costs

Accurate Across the Widest Dynamic Range

CellDrop's variable volume counting chamber will automatically keep your cells in the optimal range for the most accurate results and eliminate need to dilute or concentrate many samples.

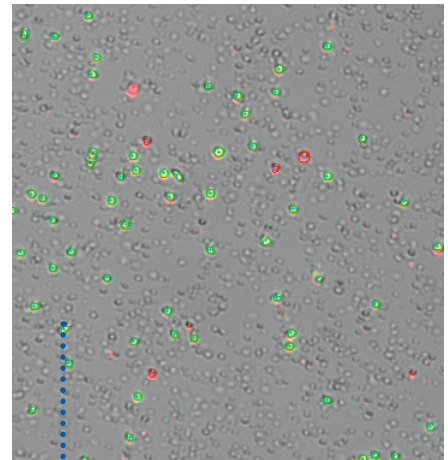

PBMCs: discriminate between live/dead and non-nucleated cells.

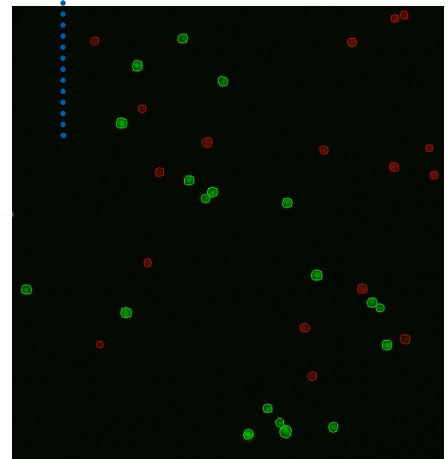

Accurate High Definition Edge Detection with regular and irregular cells, even at high density.

Reduce Plastic Use in Your Lab

An estimated 5.5 million metric tons of single-use plastic waste is generated by life science labs annually*. As plastic pollution becomes ever more pervasive, innovation to assist researchers to reduce this problem is increasingly urgent.

Accelerate your research and improve the sustainability of your lab by removing the environmental impact of disposable plastic slides.

CellDrop™




DirectPipette Technology


Drawing on the SmartPath® Technology from the award winning DS-11 Series, CellDrop brings the familiar Load, Measure & Wipe Clean functionality of microvolume spectrophotometers to cell counting. Simply lower the arm, pipette 10µL of cells and press count. Wipe clean with a dry lab wipe and the CellDrop is ready for the next sample with no concerns about carryover.

The Most Powerful, Easy-To-Use Cell Counter Platform

- Three-Channel Image Overlay
- Cell Gating and Cell Size Histograms
- Advanced Declustering Algorithms
- Unlimited Custom User Methods and Protocols
- Powerful Data Reporting
- Password-Protected User Accounts

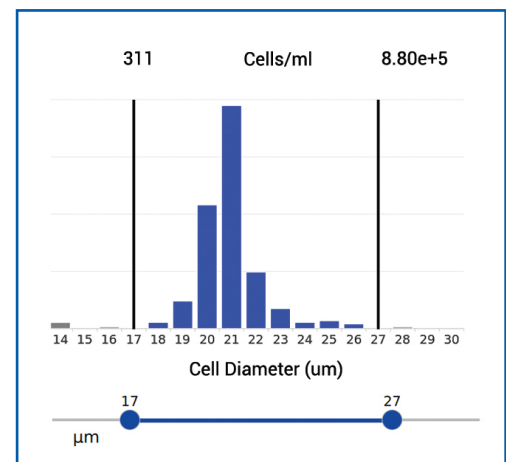

Cell Size Histograms enable simple size gating and exclusion other cell populations or debris

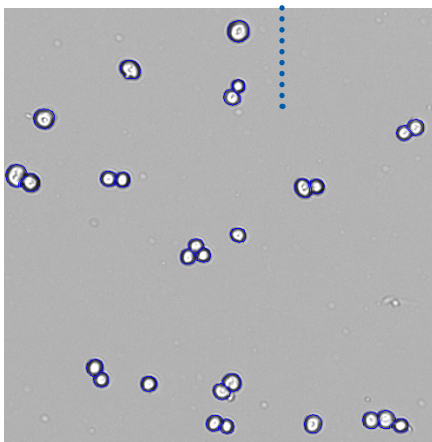

Advanced Declustering Algorithms embedded in EasyApps software improve accuracy and reproducibility even with the most challenging cells

● Dual Fluorescence and Brightfield Options

The CellDrop offers the ultimate in flexibility to meet the needs of today's leading cell culture labs. The CellDrop FL combines dual fluorescence optics with brightfield illumination to provide the widest feature set for current and future application needs. The CellDrop family also includes the CellDrop BF, Brightfield Cell Counter, ideal for labs to automate and reduce the cost of trypan blue counting.

* Mauricio A. Urbina, Andrew J. R. Watts, Erin E. Reardon. Environment: Labs should cut plastic waste too. Nature, 2015; 528 (7583): 479 DOI: 10.1038/528479c

Hazardous Samples?

The small footprint and onboard processing allow the CellDrop to fit inside most flow hoods. For applications that require the cells to be contained when reading, the CellDrop family has the option to use common disposable plastic slides and reusable slides without the need for additional adaptors.



21 CFR Part 11 Compliance Ready

EasyApps Secure provides an optional suite of software controls that allow regulated GxP facilities to easily add CellDrop to their cell counting workflow. The software is fully integrated within the on-board operating system and includes a range of features essential to ensuring compliance, including:



- Password-Protected System Access
- Integrated Electronic Signature Controls
- Comprehensive User Account Management
- Secure Audit Trail Recording
- Advanced Data Handling and Export Tools

Other Products from DeNovix



DS-11 Series

1 μ L Spectrophotometer / Fluorometer



QFX

4-Channel Fluorometer



DS-C

Cuvette UV-Vis Spectrophotometer

Specifications

Dynamic Range	7 x 10 ² to 2.5 x 10 ⁷ cells/ml	Images	FL: 2048 x 1536 px (3.15 MP) BF: 2592 x 1944 px (5 MP) with overlay capabilities
Sample Volume	5.0 μ L (high density), 10 μ L (standard), 40 μ L (low density)	Focus	Software controlled course and fine adjustment
Measurement Speed	Brightfield 3 seconds AO/PI - 8.5 seconds	Detector	FL: Sony IMX265 BF: ON Semiconductor MT9P031
CellDrop FL Modes	Dual Fluorescence Single Fluorescence Fluorescence + Brightfield	Processor	NVIDIA Tegra K1 with Quadcore ARM Cortex-A15
CellDrop BF Modes	Brightfield only	Connectivity	Wi-Fi, Ethernet, HDMI, 3 USB ports
Sample Surfaces	Optical Sapphire	Footprint (LxWxH)	37cm X 21 cm X 18 cm
Brightfield Illumination	LED 530 nm	Weight	8 kg
Fluorescence Illumination	LED 470 nm	Operating Voltage	12 VDC
Emission Filters	AO 525 nm +/- 25 nm, PI 645 nm +/- 37 nm	Approvals	UL/CSA, CE, FCC, Japan CAB
Gesture Recognition	Multipoint touch, swipe, pinch	Manufacture Location	USA
Display	High Definition 1280 x 800 7" color display	Warranty	2 Years
Glove Compatibility	All common lab gloves	Colors	Arctic White
Camera	FL: Basler acA2040-55um; BF: Basler daA2500-14um	Internal Storage	120 GB solid state drive - upgradable to 1 TB
		Accessories	Barcode reader, keyboard, mouse